

Manual Patient Lift 71910



Weight Capacity
450lbs / 204kg



Hydraulic pump
(thrust up to 9000N)

Comes standard with 5" casters
(rear wheel with lock)

Product Description


Patient lift is an assistive device, suitable for transfer people such as Wheelchair Users, Bariatric Patient, Bedridden patients, Disabled, Handicap, lower limb fracture patients, elderly people with poor ability to walk to a wheelchair.

Transfers may be to and from beds, chairs, floor to bed, lateral transfers, bathing, and toileting. The equipment is compact in design and versatile, suitable for moving all types of patients. Lightweight and easy to disassemble, convenient for fast transportation and installation. Can make patients and caregivers feel safe, comfortable and stable.

Features:

1. Manual Hydraulic Patient Lift
2. Operation optimizes lift capabilities to require less physical effort for raising patient
3. High-performance hydraulics for safe and gradual movement
4. Sturdy, silver vein steel construction and 6-point cradle design provides maximum strength to raise and lower individuals weighing up to 450 lbs.
5. Standard, 3" or 5" casters are offset for greater stability and smooth operation
6. Easy-to-operate caster breaks provide additional safety and security
7. 6-Point Swivel Bar can accommodate 2 or 4 sling straps plus 2 chain connections

Applications:

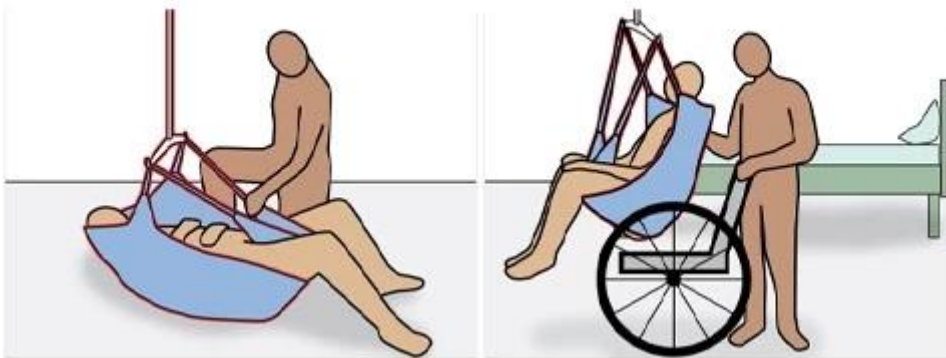



Lifting from floor
Toileting/Bathing

Bed-Chair transfer

Detailed Description



Secure and Safe:

Six-point swivel bar attachment easily adapts
to all sling styles and positions
(can be used with 4 or 6 point sling).



6.5 inch standard base design, the base can be installed under bed, wheelchair or furniture. Expand base of lifter to widest possible position to maximize stability.



The Deluxe Hydraulic Patient Lift can be used with chains or slings with straps. The four-point swivel bar attachment easily adapts to most sling styles and positions.

Product Comparison:



	Hydraulic/Manual Patient Lift 71910	Battery Powered Patient Lift 71920	Bariatric Battery Powered Lift 71922	Electric Patient Lift71940
 Weight Capacity	450 lbs/204kg	450 lbs/204kg	600 lbs/272kg	400 lbs/181kg
 Net Weight	32.7kg	49kg	59.9kg	39.9kg
 Gross weight	38.1kg	55.3kg	67.1kg	48.5kg
 Base Height	6.5"	4.75"	4.5"	4.4"
 Boom Height	17.5"-72"	27.5"-76"	32" - 75.5"	20" - 65"
 Casters	5"	3" Front/5" Rear	3" Front/4" Rear	4"
 Min and Max Base Width	42" / 22.5"	41" / 24"	45.75" / 25"	33" / 23"
 Weight Scale	/	Optional	Optional	/