

# SPECIFICATIONS

|                    |                            |
|--------------------|----------------------------|
| Boom Height        | 17.5" - 72"                |
| Base Width         | 42"(open) / 22.5" (closed) |
| Base Height        | 6.5"                       |
| Casters            | 5"                         |
| Net Weight         | 72 lbs/32.7kg              |
| Gross Weight       | 84 lbs/38.1kg              |
| Weight Capacity    | 450 lbs/204kg              |
| Package Dimensions | 1250*630*342mm             |



**450 lbs**

weight capacity

**APPLICATIONS**



## DETAILED DESCRIPTION

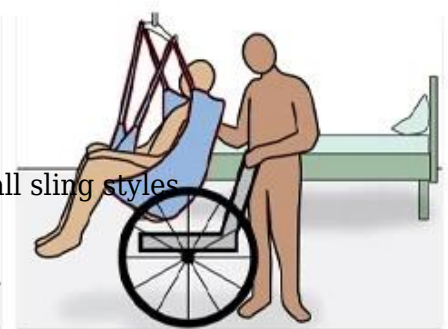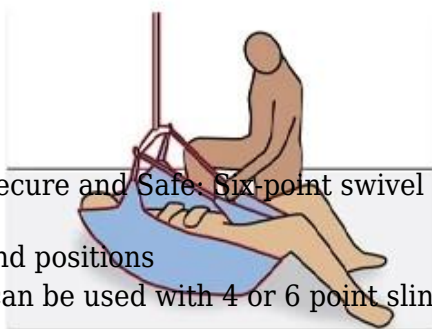

Secure and Safe: Six-point swivel bar attachment easily adapts to all sling styles and positions (can be used with 4 or 6 point sling).

You will also want to give some thought as to the best sling for the patient, in terms of support and security.



The lift utilizes high-performance hydraulics to ensure safe and gradual movement.

The maximum thrust is 9000N.

Pump handle can rotate from side to side for convenience of Attendant.

Swivel bar and pump handle are padded for comfort.



6.5 inch standard base design, the base can be installed under bed, wheelchair or furniture.

Expand base of lifter to widest possible position to maximize stability.



The 71910 offers a lift range of 17.5" to 72", so not only can it lift patients from the floor, but it can also facilitate transfers from other surfaces, like chairs and hospital beds.

## **SIMILAR PRODUCTS**

